



Empfohlene

Reinraumklassen

ISO 3|4|5|6|7|8|9

GMP C|D



Halyard Nitril-Handschuhe PUREZERO HG3 blau #CLN9031SM

pure¹¹-Nr.: 1105010, Marke: Halyard

Eigenschaften

- Marke: Halyard
- Handschuhtyp: Dünnsfilm
- Länge in cm: 30,5 cm
- Chemikalienbeständigkeit - Typ: Typ C
- Puderfrei
- Material: Nitril
- Texturierte Fingerspitzen
- Zytostatikageeignet
- Lebensmittelkonformität
- Silikonfrei
- Vulkanisationsbeschleunigerfrei
- Verpackungsform: Beutel
- Antistatisch
- AQL (Acceptable Quality Level)-Wert: 0,65
- Chemikalienbeständig gegen Isopropanol (70%): Level 1 (10-30min)
- Länge in Inches: 12 In
- Latexfrei
- Reißfestigkeit EN 455-2 ASTM in MPa: 20-30
- Schutz vor Blut und Körpersekreten ISO 16604:2004

- Passform Hand: beidhändig
- Viren-/Mikroorganismenschutz EN ISO 374-5:2016
- Wandstärke Mittelfinger in mm: 0,1 mm

Material

- Nitril

Verpackung

- 1500STK

Produktvarianten

pure¹¹-Nr.: 1105010BLS, Halyard Nitril-Handschuhe PUREZERO HG3 blau #CLN9031SM

Farbe: Blau; Größe: S / VE: 1500STK

pure¹¹-Nr.: 1105010BLL, Halyard Nitril-Handschuhe PUREZERO HG3 blau #CLN9031LG

Farbe: Blau; Größe: L / VE: 1500STK

pure¹¹-Nr.: 1105010BLM, Halyard Nitril-Handschuhe PUREZERO HG3 blau #CLN9031MD

Farbe: Blau; Größe: M / VE: 1500STK

pure¹¹-Nr.: 1105010BLXL, Halyard Nitril-Handschuhe PUREZERO HG3 blau #CLN9031XL


Farbe: Blau; Größe: XL / VE: 1500STK

pure¹¹-Nr.: 1105010BLXS, Halyard Nitril-Handschuhe PUREZERO HG3 blau #CLN9031XS

Farbe: Blau; Größe: XS / VE: 1500STK

GUIDE TO DONNING AND DOFFING OF NON-STERILE CLEANROOM GLOVES

DONNING GLOVES

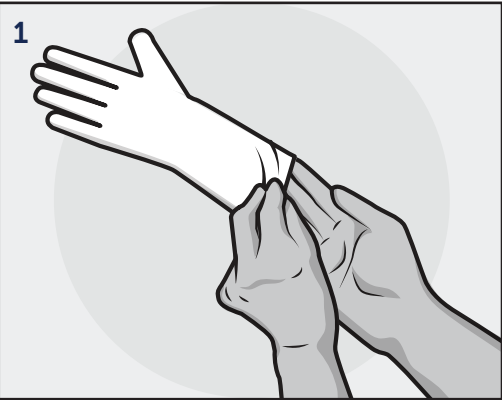


Grasp the beaded cuff of one glove and carefully pull over the opposite hand.




Gently ease the glove up over the hand and wrist.



Grasp the second glove at the edge of the beaded cuff.

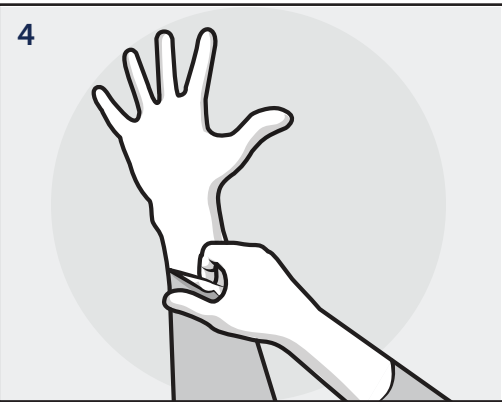


Gently ease the glove up over the hand and wrist.

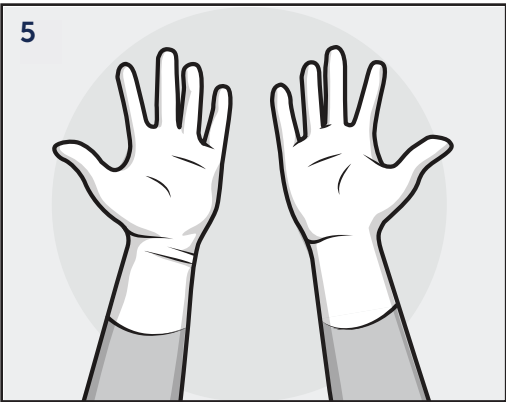


You are now ready to perform your task.


DOFFING GLOVES

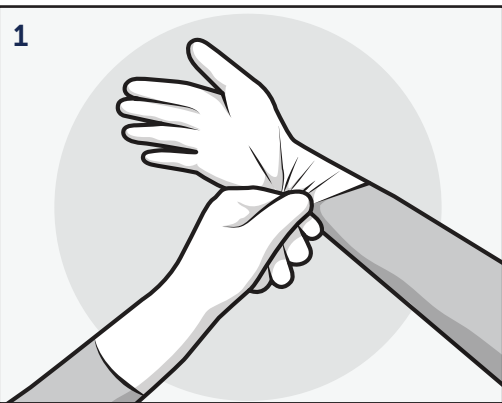

Grasp the exterior of one glove with your other gloved hand.

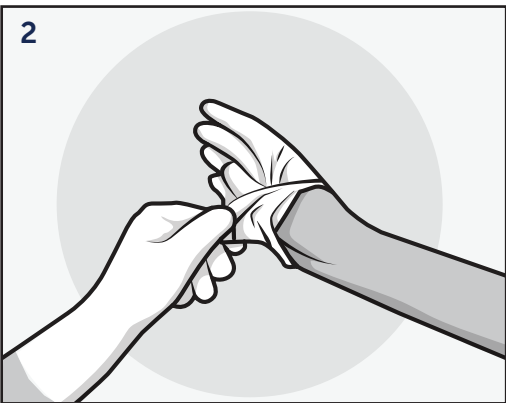


Carefully pull the glove off your hand, turning it inside-out.

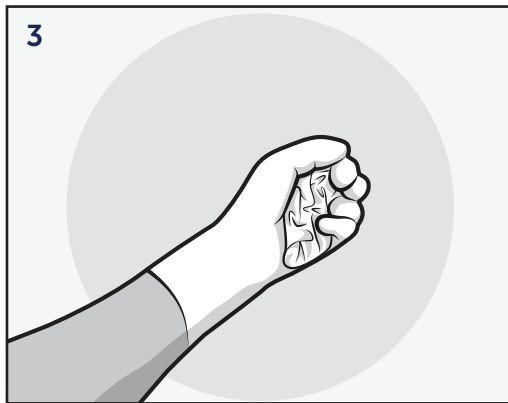


Ball the glove up and hold in your other gloved hand.



Slide your ungloved fingers into the opening of the other glove.

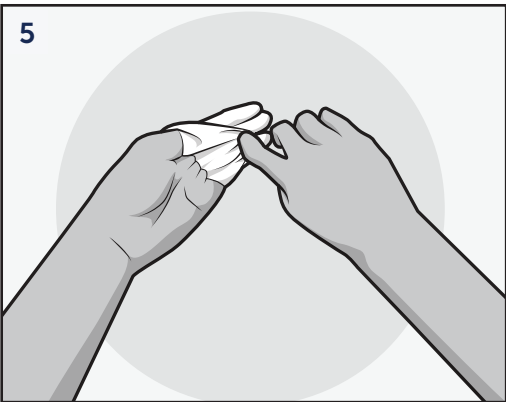

Carefully pull the glove off your hand, turning it inside-out and wrapping it around the doffed glove in your palm.

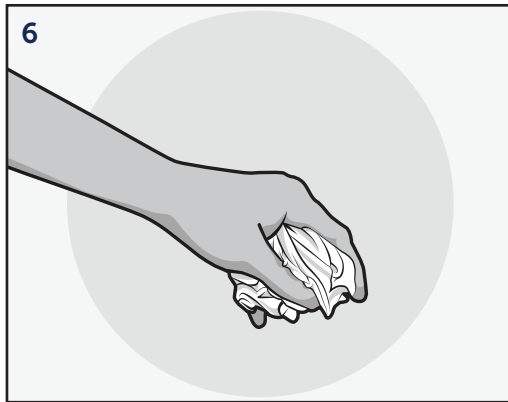

Discard appropriately.

For more information, contact your HALYARD* representative or distributor, or visit www.purezerogloves.com