


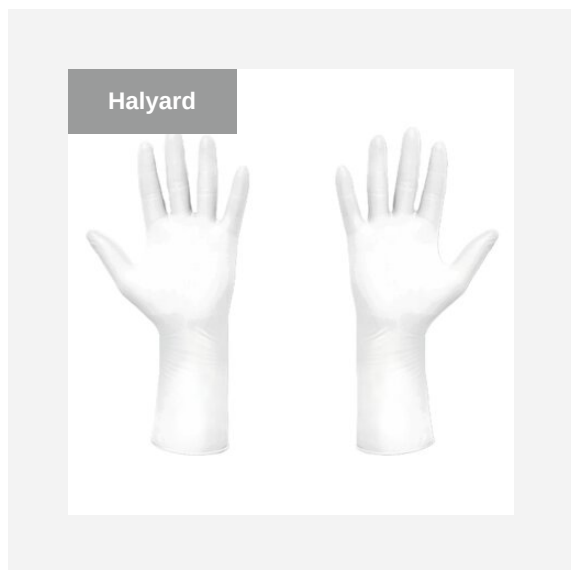

Empfohlene

Reinraumklassen

ISO 3|4|5|6|7|8|9

GMP C|D



Halyard Nitril-Handschuhe PUREZERO HG3 weiß #CLN3031SM

pure¹¹-Nr.: 1105011, Marke: Halyard

Eigenschaften

- AQL (Acceptable Quality Level)-Wert: 0,65
- Chemikalienbeständig gegen Isopropanol (100%): Level 3 (60-120min)
- Chemikalienbeständig gegen Isopropanol (70%): Level 3 (60-120min)
- Chemikalienbeständigkeit - Typ: Typ B
- Handschuhtyp: Dünnschicht
- Länge in cm: 30,5 cm
- Länge in Inches: 12 In
- Latexfrei
- Lebensmittelkonformität
- Material: Nitril
- Materialzusammensetzung: Reinmaterial
- Oberflächenbeschaffenheit: griffig
- Puderfrei
- Reißfestigkeit EN 455-2 ASTM in MPa: 20-30
- Schutz vor Blut und Körpersekreten ISO 16604:2004
- Silikonfrei
- Texturierte Fingerspitzen
- Passform Hand: beidhändig
- Verpackungsform: Beutel

- Viren-/Mikroorganismenschutz EN ISO 374-5:2016
- Vulkanisationsbeschleunigerfrei
- Wandstärke Mittelfinger in mm: 0,16 mm
- Zytostatikageeignet
- Marke: Halyard

Material

- Nitril

Verpackung

- 1000STK

Produktvarianten

pure¹¹-Nr.: 1105011WHS, Halyard Nitril-Handschuhe PUREZERO HG3 weiß #CLN3031SM

Farbe: Weiß; Größe: S / VE: 1000STK

pure¹¹-Nr.: 1105011WHL, Halyard Nitril-Handschuhe PUREZERO HG3 weiß #CLN3031LG

Farbe: Weiß; Größe: L / VE: 1000STK

pure¹¹-Nr.: 1105011WHM, Halyard Nitril-Handschuhe PUREZERO HG3 weiß #CLN3031MD

Farbe: Weiß; Größe: M / VE: 1000STK

pure¹¹-Nr.: 1105011WHXL, Halyard Nitril-Handschuhe PUREZERO HG3 weiß #CLN3031XL


Farbe: Weiß; Größe: XL / VE: 1000STK

pure¹¹-Nr.: 1105011WHXS, Halyard Nitril-Handschuhe PUREZERO HG3 weiß #CLN3031XS

Farbe: Weiß; Größe: XS / VE: 1000STK

GUIDE TO DONNING AND DOFFING OF NON-STERILE CLEANROOM GLOVES

DONNING GLOVES

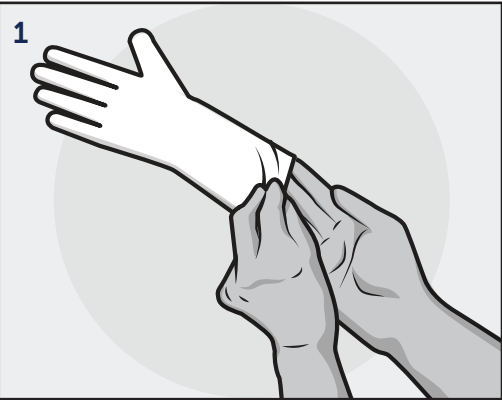


Grasp the beaded cuff of one glove and carefully pull over the opposite hand.




Gently ease the glove up over the hand and wrist.



Grasp the second glove at the edge of the beaded cuff.

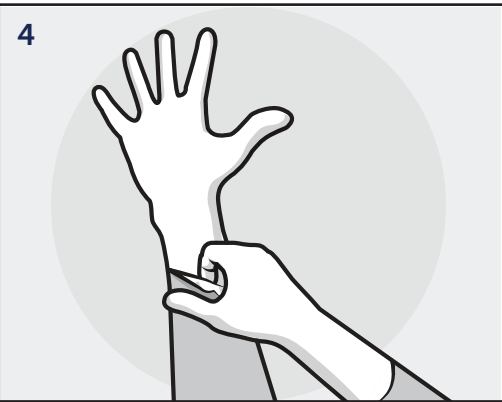


Gently ease the glove up over the hand and wrist.

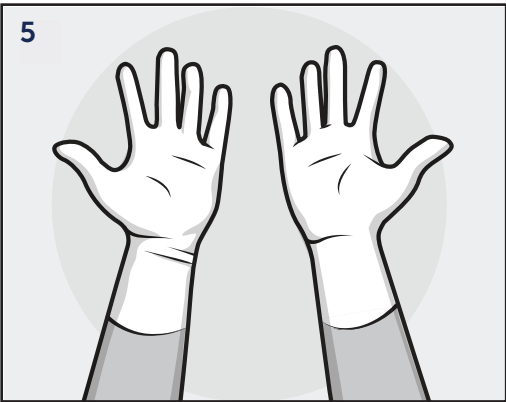


You are now ready to perform your task.


DOFFING GLOVES

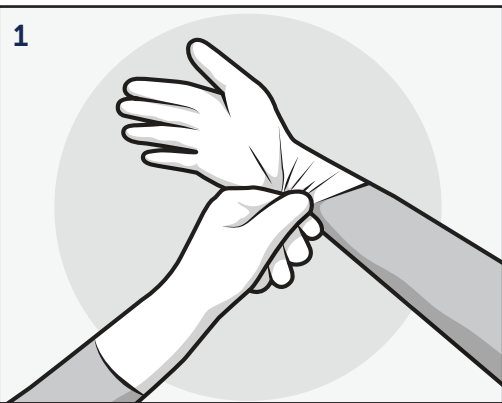

Grasp the exterior of one glove with your other gloved hand.

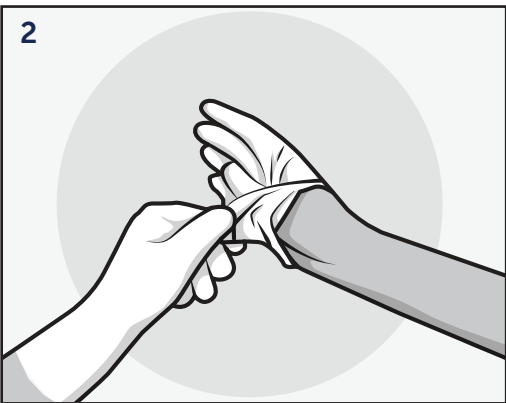


Carefully pull the glove off your hand, turning it inside-out.

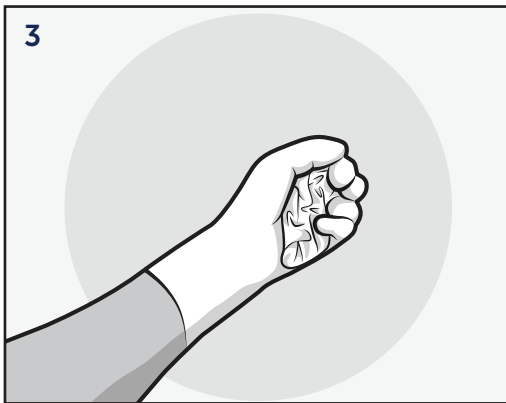


Ball the glove up and hold in your other gloved hand.



Slide your ungloved fingers into the opening of the other glove.

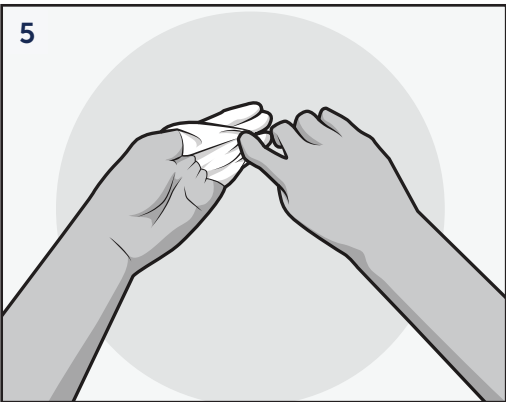

Carefully pull the glove off your hand, turning it inside-out and wrapping it around the doffed glove in your palm.

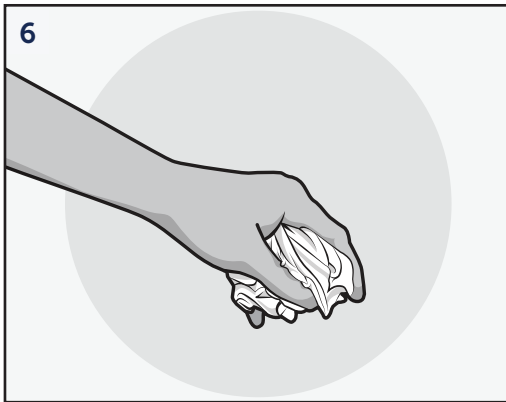

Discard appropriately.

For more information, contact your HALYARD* representative or distributor, or visit www.purezerogloves.com